



# ***LAB Online Exhibition***



## Overview

This Quick Start guide is intended to help you with the basic setup and operation of your Brasonic Ultrasonic Bath.

Consult your Operator's manual for detailed information on operation, applications, solutions, maintenance, and troubleshooting.


### WARNING

- Do not use alcohol, gasoline or flammable solutions. Doing so could cause a fire or explosion. Use only water-based solutions.
- Do not operate the unit without proper grounding.
- Do not disassemble your unit—high voltage inside the unit is dangerous.
- Do not immerse the unit in water.
- Do not operate the unit without liquid in the tank. This will cause the unit to fail.
- Do not allow the solution to drop more than 3/8 inch (1 cm) below the operating level line with the heater on. Doing so will cause the heater to fail.
- Do not use mineral acids or chlorine solutions such as chlorine bleach. These could damage the tank.
- Do not place parts or containers directly on the bottom of the tank. Use a tray or wire to suspend items. Direct placement can cause the unit to fail.

## Installation

For models **3800**, **5800**, and **8800**:


1. Remove thread protecting cap and hand tighten the drain valve.
2. Finish tightening the valve using a 21 mm or an adjustable wrench. Tighten no more than 1 full turn until the valve handle is on top. **DO NOT OVERTIGHTEN.**
3. Hand tighten the hose adaptor. Slide the drain tube over adaptor end.

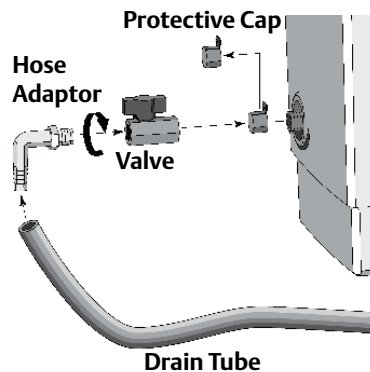

For **all** Models:

1. Plug detachable line cord into the back of the unit
2. Connect the unit to a grounded power source. (Consult data-tag on back of unit for electrical ratings.)

## M/MH Series Operation

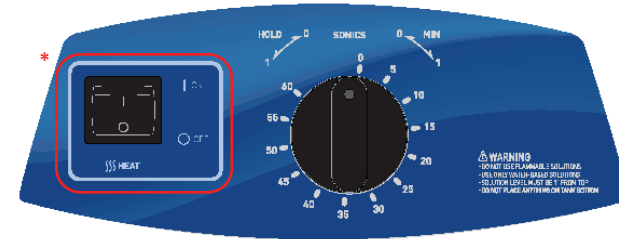

\* Heater available on MH models only


Step	Action
1	Fill tank with water close to the operating level line and add solution per the manufacturer's recommendation.
2	Turn Main Power switch on (located next to the power cord receptacle)
3	Turn the Timer Knob clockwise to 5–10 and let the unit run to allow the solution to "degas."
4	Turn HEAT switch on, if available.
5	Turn the Timer Knob clockwise to set the amount of time (0–60 minutes) you wish the items to be cleaned. Turn the Timer Knob counterclockwise to the HOLD position for continuous operation.
6	Place the items into a basket, perforated tray, or beakers in a positioning cover.
7	If using beakers or a solid tray, add cleaning solution to beakers or tray to cover the items.
8	<b>Slowly</b> lower the tray or beakers into the tank. Do not allow items to contact the tank bottom. Do not stir the solution.
9	When items are clean, <b>slowly</b> remove them from the tank.
10	Rinse the clean items with clean water and dry them, if necessary.

## CPX Series Operation



Step	Action
1	Fill tank with water close to the operating level line and add solution per the manufacturer's recommendation.
2	Turn Main Power switch on (located next to the power cord receptacle)
3	Press the On/Standby key (⏻) to turn on the unit.
4	Press Degas key (🧼) once to start the degas process. Default degas time is 5 minutes.
5	Set the amount of time you wish the items to be cleaned, or select Constant Sonics Mode: <ul style="list-style-type: none"> <li>Use up/down keys to increase/decrease cycle time (hold for quick increments/decrements).</li> <li>Pressing up key from 99 minutes or down from 1 minutes takes you to Constant Sonics Mode (∞ and "- -" on the display).</li> </ul>
6	Press the Sonics key (🔊) to activate ultrasonics.
7	Place the items into a basket, perforated tray, or beakers in a positioning cover.
8	If using beakers or a solid tray, add solution to beakers or tray to cover the items.
9	<b>Slowly</b> lower the tray or beakers into the tank. Do not allow items to contact the tank bottom. Do not stir the solution.
10	When items are clean, <b>slowly</b> remove them from the tank.
11	Rinse clean items with clean, warm water and dry, if necessary.

## CPXH Series Operation

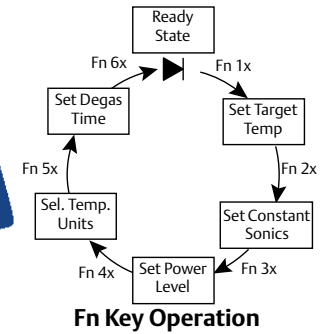

Step	Action
1	Fill tank with water close to the operating level line and add solution per the manufacturer's recommendation.
2	Turn Main Power switch on (located next to the power cord receptacle)
3	Press the On/Standby key (⏻) to turn on the unit.
4	Press Degas key (🧼) once to start the degas process. Default degas time is 5 minutes.
5	Set the amount of time you wish the items to be cleaned: <ul style="list-style-type: none"> <li>Use up/down keys to increase/decrease cycle time (hold for quick increments/decrements).</li> </ul>
6	Set the tank temperature: <ul style="list-style-type: none"> <li>Press the Fn key until the Set Temperature icon (🌡️) flashes.</li> <li>Press the up/down keys to alter the setting to the tank temperature you wish to maintain. While blinking, press the Fn key 5x or wait 15 seconds to return to the ready state.</li> <li>Press the HEAT key (🔥) once to activate heat. The Heat icon (🔥) appears.</li> </ul>
7	Press the Sonics key (🔊) to activate ultrasonics.
8	Place the items into a basket, perforated tray, or beakers in a positioning cover.
9	If using beakers or a solid tray, add cleaning solution to beakers or tray to cover the items.
10	<b>Slowly</b> lower the tray or beakers into the tank. Do not allow items to contact the tank bottom. Do not stir the solution.
11	When items are clean, <b>slowly</b> remove them from the tank.
12	Rinse clean items with clean, warm water and dry, if necessary.