

Ultracompact water activity kit.

- Low-maintenance sensor and readout in a single unit.
- Perform quick checks, measure a_w in 5 minutes or less.
- Easy-to-read display shows water activity and sample temperature together.
- Long 3-year battery life. No adapter or external power supply required.
- Remarkably portable, and lightweight.
- One year warranty—parts and labor.
Extended warranty available.
- Easy calibration saves time.
- Included: protective case, cups, solutions, quick-start guide, and manual.



■ Designed for quick checks by health inspectors, food scientists, and QC managers in foods, pharmaceuticals, and cosmetics.

Water activity—four inches sleek and four ounces slim.

**Pawkit
AquaLab**

When you can't bring the lab with you. The new portable Pawkit water activity kit is a lightweight, space-saving device for any on-the-go Quality Control Professional.

Four inches.

This sleek device is compact, lightweight, portable and discreet. At four inches long, it takes up no more room than a small calculator.

Four ounces.

The ultra-compact Pawkit delivers $\pm 0.02 a_w$ measurements in 5 minutes, perfect for quick checks at facilities without laboratories. ■

DECAGON

Post Office Box 835
Pullman, Washington 99163

**Call today.
800-755-2751**

fax 509-332-5158
pawkit@decagon.com
www.decagon.com/pawkit/



SPECIFICATIONS

Pawkit

Operation

environment:

5 to 50°C (41 to 122°F).

0 to 90% relative humidity
(non-condensing).

Sensor:

Dielectric humidity sensor.

Measurement speed:

5 minutes.

Power:

2 Li Coin cells (3V).

Battery life:

3 years typical.

Display:

6-digit LCD display with
symbols.

3-digit a_w , 3-digit
temperature.

Accuracy:

$\pm 0.02 a_w$.

Range:

0 to 1.0 a_w .

Resolution:

$\pm 0.01 a_w$.

Weight:

115g (4oz).

Case dimensions:

3.5 x 4 inches (oval).

Warranty:

One year-parts and labor.
Extended warranty
available.

Case:

Stainless steel with
silicone elastomer base.

Sample dish

capacity:

7ml recommended
(15ml full).

Calibration

standards:

6.0m NaCl - 0.76 a_w .

13.41m LiCl - 0.25 a_w .


Carrying case:

Clamp-hinge rigid shell
with foam insert and
external vinyl protection.
(9.75" x 6" x 2").



■ When your
work takes
you away
from the lab,
take along
Pawkit.

*Measuring water activity
away from the lab has
never been so easy.*

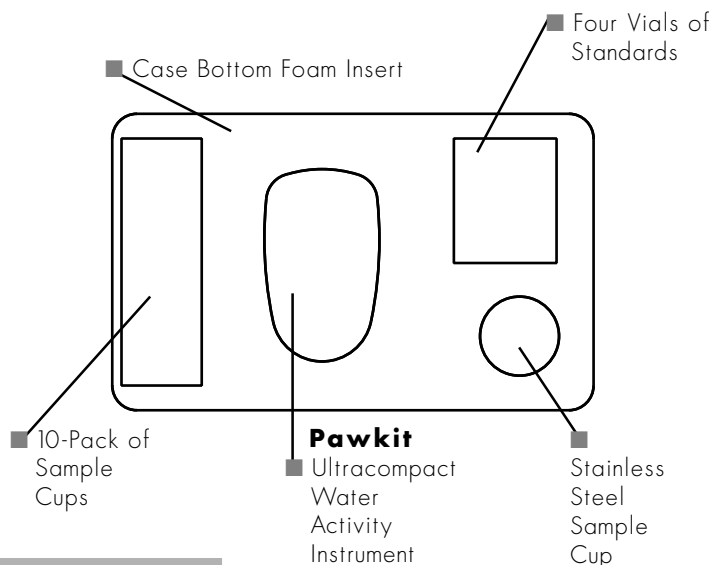

Additional case kit

accessories (not shown):

- Extra 100 disposable sample cups.
- Quick-start guide & user's manual.

DECAGON

Post Office Box 835
Pullman, Washington 99163

800-755-2751

fax 509-332-5158

pawkit@decagon.com

www.decagon.com/pawkit/

DS14PA-10

©2000 Decagon Devices, Inc. Printed in USA.