


Our most compact portable planetary centrifugal mixer

AR-100

AR-100 THINKY MIXER

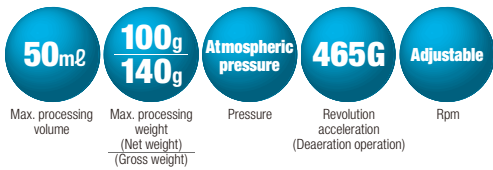

The AR-100, the smallest model of the THINKY Mixer Series, is designed for user-friendliness and the repeated processing of small batches of material. The space-saving and compact design is ideal for basic experiments conducted by researchers and engineers. By adopting disposable containers the AR-100 is perfect for experiments and research requiring frequent use of the unit.

Features

- About 460 Gs of acceleration generated by rotation and revolution speed allows powerful simultaneous mixing and deaeration.
- By changing the rpm, optimal mixing of any material is ensured.
- Timer setting with 5 recipe profiles.
- Easy-to-open and close sliding lid.
- A stroboscope allows observation of the material in a stationary state.

Unit dimensions:

H328 × W250 × D250 (mm)

Unit weight: Approx. 15 kg

Max. processing volume with a 100 ml disposable container:

50 ml / 100 g (net)

50 ml / 140 g (gross)



● Container holder, Balance adjustment dial



● Control panel



● Accessories: Containers

50ml

100g
140g

Atmospheric pressure

465G

Adjustable