


粉體比重試驗器 Powder Tapping Tester

◎量筒無附屬

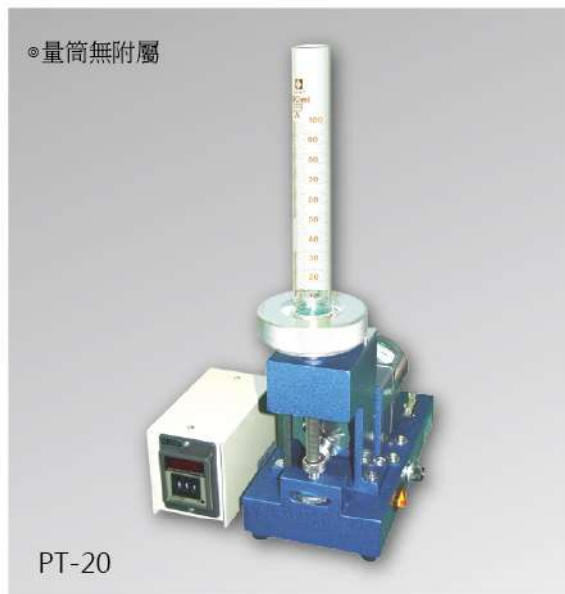

PT-20

特點：

本機是藉由規律地拍打含有粉末的測量桶來達成，當使用者選擇充填的量後，便可開始拍擊的動作。當產生改變時，讀數便會產生。此種拍擊是藉由測量筒的上下振動，在特定的重量及指定的距離來完成。在拍擊的過程中，會優先地將任何可分開的團塊物，在向下拍擊的同時將之最小化。而使用透明的塑膠填充管則有助於觀察及讀出測量值。

Feature :

This machine is to achieve the effect by a regularly tapping on the powder filled cylinder. When the user determines the padded volume, the tapping function can be started. As the density of the powder is changed, the readings will be followed by changes. Such a tapping function is completed by the vertical vibration of the cylinder with a specific weight and a specified distance. Any nodules can be first separated into the smallest particles, and the transparent plastic tube is very convenient to observe and read out the measured values.

Model	PT-20
振幅 Stroke (mm)	10 ~ 45 Adjustable
次數 Tapping Times (min)	0 ~ 100
容量 Tapping Tube Capacity (ml)	10 ~ 250
量筒座直徑 Tapping Base Diameter (mm)	< Ø70
數位顯示次數 Digital Times	1 ~ 999
馬達 Motor	AC 25W
外型尺寸 Outer Size (W x L x Hmm)	220 x 270 x 240
重量 Weight (kg)	6.5
電源 Power	AC 220V, 50/60Hz