

Simple, Robust and Mobile Pro2Go pH Meter



Versatile and User-Friendly

Pro2Go™ offers easy set-up from its intuitive menu (in 10 languages) so that users can simply take it out of the box and operate it with no need to study complicated manuals. The meter's sleek, ergonomic design makes it easy to handle.




Durable for Industrial Use

A pH meter used in the field needs to stand up to your industrial environment. Pro2Go provides IP67 ingress protection from dust and water damage. It has been drop tested to ensure that it can withstand harsh and demanding applications.



Easy Data Capture and Management

Pro2Go collects data from measurements and stores it in its internal memory. It can store up to 2000 data sets and the pH meter can connect via micro-USB to a computer for data transfer and for powering.

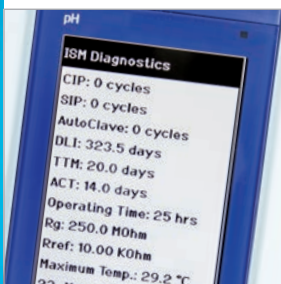

ISM Diagnostics in a Handheld Device

The Pro2Go portable pH meter is fully compatible with both analog pH sensors and pH sensors with digital Intelligent Sensor Management (ISM®) technology. This means it has full sensor diagnostics including DLI, ACT and TTM.



Convenient pH Measurement for a Wide Range of Applications

For periodic pH or redox measurements, a mobile meter is the ideal choice. Designed for laboratory and industrial applications, Pro2Go's intuitive menu ensures out-of-the-box operation for anyone.

The meter's ergonomic design allows one-handed operation for both big and small hands. Its light weight makes measurement easy and convenient, even over lengthy periods of repeated measurement.

Pro2Go is compatible with both analog and digital ISM sensors. Its ISM features include Plug and Measure for error-free setup, and display of sensor diagnostics including DLI, ACT and TTM, so operators know when sensor calibration or replacement will be required.

Technical Data of the Pro2Go pH Meter

Measurement parameters	pH, mV and temperature
Supply voltage batteries	4 × LR6/AA 1.5V alkaline or 4 × HR6/AA 1.2V NiMH rechargeable
Sensor type	Analog and ISM sensors
Supply voltage (USB powered)	Connection: Micro-USB Rating: 5VDC, 100mA
pH measuring range	-2 to 20 pH
User interface	Graphic LCD display
pH resolution	0.001 / 0.01 / 0.1 (selectable)
PC connection	Micro-USB for data transfer and power
Battery life (standby)	200 to 250 hours
Ambient temperature	0 ... 40 °C (32 to 104 °F)
Weight	290 g (0.64 lb)
Data logger	2000 datasets (GLP conform)
Enclosure rating	IP 67

Ordering information

Pro2Go portable pH Meter	Order No.
Pro2Go portable pH Meter including USB cable, sensor cable AK9-BNC/RCA for ISM sensors, rubber holster, wrist strap, CD with documentation and software, Declaration of conformity, Test certificate	30 386 271

Accessoires	Order No.
Rubber holster	30 487 344
USB cable for PC connection	30 487 345
Power adapter for USB cable (to operate instrument without batteries)	30 487 346
Sensor cable AK9-BNC/RCA for ISM sensors	30 487 466
EasyDirect pH PC software	free download

ISM is a registered trademark of the METTLER TOLEDO Group in Switzerland, Brazil, USA, China, European Union, South Korea, Russia and Singapore. Pro2Go is a trademark of the METTLER TOLEDO Group.