

Rotary Evaporator

The SCIOGEX SCI100-S Digital Rotary Evaporator is an excellent choice for evaporation, distillation and separation of chemicals. It can be used as a multifunctional system in combination with vacuum pump, vacuum controller and circulating refrigerators to provide excellent distillation solutions in many laboratory applications.

Features:

- 5L heating bath with a wide temperature range (room temp. to 180°C), independent temperature control, make it can be used separately. Water/oil heating mode can be changed only through a switch.
- Manual lift combined with auxiliary lift for precise positioning of the glassware.
- PID control ensures high temperature accuracy at $\pm 1^\circ\text{C}$ (water).
- Overheating protection temperature at 220°C
- Boil-dry protection, automatically power off if heating without water/oil in the heating bath.
- Speed range from 20 to 200rpm, timing interval operation in clockwise and anticlockwise directions for drying process.
- Patented condenser (cooling surface 1200/1700cm² for selection), double helix condensing tube + center circular arc design to speed up the flow rate of liquid with excellent cooling effect.
- Ejection mechanism ensures easy exchange of evaporating flask
- Adjustable immersion angle
- Patented double spring sealing ring which is made of PTFE provides an excellent sealing performance
- Remote function provides PC control and data transmission.
- Compatible with the entire range of RE100-Pro glassware.



USB connector



Chemical resistance




Overheating protection



Compatible with all range of glassware

SCI100-S LED Digital Rotary Evaporator

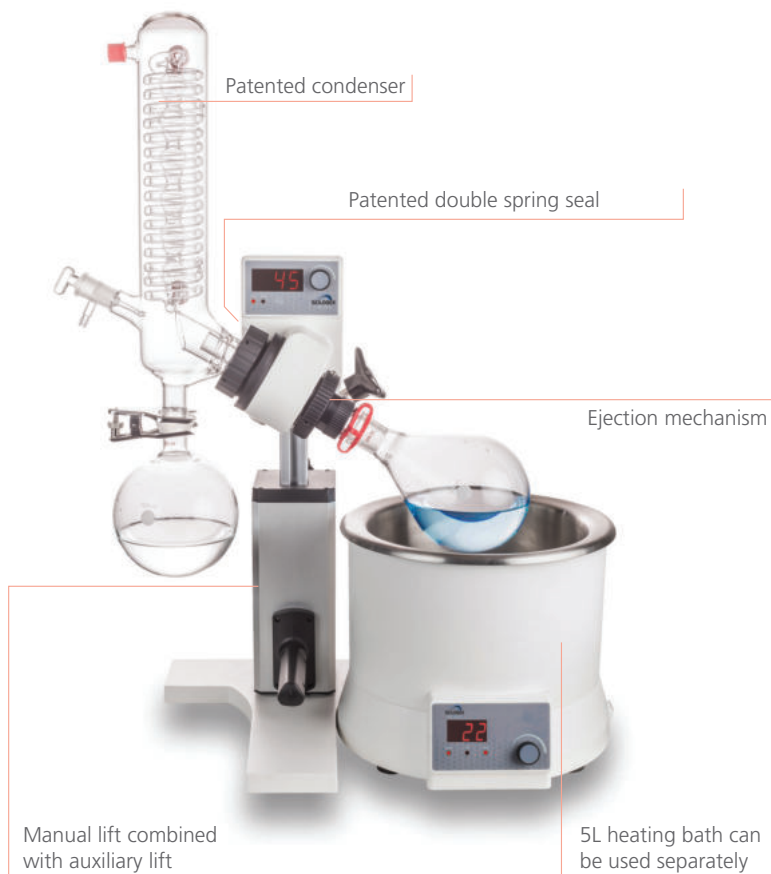

Specifications

	SCI100-S
Motor Type	Brushless DC motor
Speed Range	20-200rpm
Display	LED (speed, temperature, time)
Clockwise and Anti-clockwise	Yes
Heating Temperature Range	Room temp. to 180°C
Control Accuracy	water: $\pm 1^\circ\text{C}$ oil: $\pm 3^\circ\text{C}$
Heating Power	1010 W
Stroke Displacement	Manual 110mm + auxiliary 100mm
Interval Time Setting Range	1~999s
Dimension [D × W × H]	Main Unit: 440 × 320 × 450mm Heating Bath: 300 × 300 × 240mm
Weight	Main Unit: 7kg Heating Bath: 3kg
Permissible Ambient Temperature	5~40°C
Permissible Relative Humidity	80%RH
Protection Class	IP20
USB Interface	Yes
Voltage/Frequency	100~120/200~240V 50/60 Hz
Power	1100 W

Condenser



Vertical

Cat. No.	Glassware type	Volume	Cooling surface
18900165	Evaporating flask(NS29/32) Receiving flask(KS35/20)	1000mL	1700cm ²
18900173	Evaporating flask(NS24/40) Receiving flask(KS35/20)	1000mL	1700cm ²
18900166	(COATED)Evaporatingflask(NS29/32) Receiving flask (KS35/20)	1000mL	1700cm ²
18900174	(COATED)Evaporatingflask(NS24/40) Receiving flask (KS35/20)	1000mL	1700cm ²
18202410	evaporating flask (NS 29/32) Receiving flask (KS 35/20)	1000mL	1200cm ²



Diagonal

Cat. No.	Glassware type	Volume	Cooling surface
18900163	Evaporating flask(NS29/32)Receiving flask(KS35/20)	1000mL	1700cm ²
18900171	Evaporating flask(NS24/40)Receiving flask(KS35/20)	1000mL	1700cm ²
18900164	(COATED)Evaporatingflask(NS29/32) Receiving flask (KS 35/20)	1000mL	1700cm ²
18900172	(COATED)Evaporatingflask(NS24/40) Receiving flask (KS35/20)	1000mL	1700cm ²



Heating Bath

Cat. No.	Description	Volume
18900201	Rotary Evaporator, 220V	5L
18900202	Rotary Evaporator, 110V	5L
18101888	LED heating bath, 220V	5L
18101936	LED heating bath, 110V	5L



Rotary Evaporator | Accessories

Glassware



Evaporating flask

Cat. No.	Volume	Size
18300116	50mL	NS 29/32
18300117	100mL	NS 29/32
18300118	250mL	NS 29/32
18300119	500mL	NS 29/32
18300120	1000mL	NS 29/32
18300121	2000mL	NS 29/32
18300655	50mL	NS24/40
18300201	100mL	NS 24/40
18300202	250mL	NS 24/40
18300203	500mL	NS 24/40
18300204	1000mL	NS 24/40
18300205	2000mL	NS 24/40



Receiving flask

Cat. No.	Volume	Size
18300107	100ml	KS 35/20
18300108	250ml	KS 35/20
18300109	500ml	KS 35/20
18300110	1000ml	KS 35/20
18300111	2000ml	KS 35/20



Distilling spider with 5 flasks(without flask)

Cat. No.	Volume	Size
18900194	/	NS 29/32, NS24/29
18900177	/	NS 29/42, NS24/40

Foam brake

Cat. No.	Volume	Size
18300093	250ml	NS 29/32
18300207	250ml	NS 24/40



Vapor tube

Cat. No.	Volume	Size
18100998	/	NS 29/32
18300206	/	NS 24/40



Connector

Cat. No.	Volume	Size
18300094	/	NS 29/32, NS24/29
18300095	/	NS 29/32, NS19/26
18300096	/	NS 29/32, NS14/23
18300199	/	NS24/40, NS19/22
18300200	/	NS 24/40, NS14/20
18900193	/	NS 24/40, NS 29/42