

COLOR METER

ZE 6000



Large bright and easy-to-read screen

Dialogue type simple operation

Availability for versatile samples

solid, liquid, powder, pellet, film, paste, etc.

For long time use and high accuracy measurements

Flicker photometric system



OPTION

- ◆ Harzen number color number (APHA) and Gardner color can be measured.
- ◆ Fluorescence intensity can be measured using a UV ray protective filter.
- ◆ Windows-based color management software (called ColorMate 5) can be used. Using it, you can view on the screen of your personal computer chromaticity coordinate graph, transition graph and a list of data, etc.



ISO 9001 certificate: Obtained in 2000

NIPPON DENSHOKU

COLOR METER

ZE 6000

Transmittance measurement




Measurement of liquid samples in square cells



Measurement of transparent samples by setting on a sample clamp

Reflectance measurement

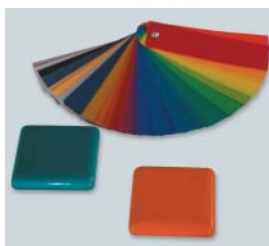

Measurement by setting color samples on a measuring stage



Measurement of powder and paste in round cells



Measurement of small samples by changing sample stages for different measuring diameters

OPTION

- ◆ Harzen color number (APHA) and Gardner color can be measured.
- ◆ Fluorescence intensity can be measured using a UV ray protective filter.
- ◆ Windows-based color management software (called ColorMate 5) can be used. Using it, you can view on the screen of your personal computer chromaticity coordinate graph, transition graph and a list of data, etc.

Features


- Both of reflectance and transmittance measurements are possible by this one unit. Its is in fact all-in-one type measuring instrument enabling color measurements of various shaped samples such as solid, liquid, powder, pellet and film, etc.
- Thanks to its flicker photometric system, measuring stability for a long time and data acquisition in high accuracy have been made possible.
- Munsell color system (HVC) can be selected from one of the C/2° and D65/2° conditions.
- Measuring diameter for reflectance can be selected from 6mm ϕ , 10mm ϕ or 30mm ϕ , with a possibility to measure even small samples. (The diameter as small as 3mm is also available on your special order, which will be most appropriate for measurement of tablet pharmaceuticals.)
- This model is equipped with a large bright and easy-to-read display.
- Color data management is available on your personal computer using the RS-232C interface.
- This model is equipped with a foot switch as one of the standard accessories, which allows you to make measurements even while you are holding a sample by both hands.



Large bright and easy-to-read screen



Display of measurement values

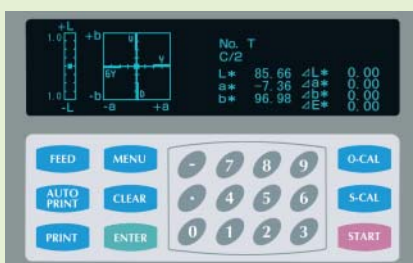

Display of color deviation criterion chart

Foot switch

Its allows you to make measurements even while you are holding a sample by both hand.



Examples of applications for different industries

Food industry

Color difference control for tomato juice, orange juice, mayonnaise, ketchup, miso, coffee, and ham

Fishing industry

Control for fish freshness and minced fish color

Cosmetics industry

Color check for foundation, lipstick, etc. and for demonstration use at shops

Dentistry

Color check for teeth, gums, and dental material

School and hospital

Educational use for housekeeping class, color control at cooking, and color check of skin

Pharmaceutical industry

Color shading control for tablets, granules, and powdered pharmaceuticals

Paper manufacturing industry

Color management for paper and pulp

Household and electric appliance industry

Color shading control for refrigerators, washing machines, cleaners, etc., and color check of TV sets, and LCD panels

Automobile and motorcycle industries

Color management for automobile and motorcycle painting

Metal industry

Color control for printing on aluminum and metal, and for general color printing

Dyeing and apparel industries

Color shading control for dyed or sewn fabric

Paint and ink industries

Color difference control for inks and pigments

Printing industry

Density analysis and color management for printed matter

Resin and plastics industries

Color management for pellets, resin sheets, moldings, and films

Petrochemical industry

Color check for petroleum

Agriculture and livestock industry

Color management for meat and processed products

Specifications

Illumination and light receiving conditions	Reflectance : Condition d (n-d) Transmittance : Condition e (n-d) Both based on JIS Z-8722
Measuring system	XYZ tristimulus direct reading
Measuring method	Double beam (flicker photometric system)
Measuring area	Reflectance : 30mm ϕ , 10mm ϕ , and 6mm ϕ (Small diameter is available in option.) Transmittance : 30mm ϕ
Light source	Halogen lamp, 12V 50W
Light receiving element	Silicon photocell quick response type
Measuring accuracy	Repeatability using a standard white plate: · ΔE^* - Standard deviation within 0.02 At the conditions of 30 measurements at intervals of 10 seconds, 15 minutes after turning on the power
Light source	Selectable from C/2, D65/2, and D65/10
Display	Number of dots : 256 $\times$ 64 dots Display area : 166.25mm $\times$ 41.45mm
Display items	L*a*b*, Δ L*a*b*, Δ E*, Lab, Δ Lab, Δ E, XYZ, Yxy, YI(E313), YI(D1925), W(CIE), W(Lab), WB, HVC, L*C*h*, Δ L*C*h*, Δ E94, Δ ECMC(1:1), Δ ECMC(2:1), LCh, Δ LCH, OD, color deviation criterion chart, Option (APHA and GARDNER) Any of the above color systems can be arbitrarily selected to display on the screen.
Printer	Thermal line dot printing system, thermal printer
Print items	Sample No., light source/field of view, data/time, L*a*b*, Δ L*a*b*, Δ E*, Lab, Δ Lab, Δ E, XYZ, Yxy, YI(E313), YI(D1925), W(CIE), W(Lab), WB, HVC, L*C*h*, Δ L*C*h*, Δ E94, Δ ECMC(1:1), Δ ECMC(2:1), LCh, Δ LCH, OD, color deviation criterion chart, Option (APHA and GARDNER) Any of the above items can be arbitrarily selected for printing.
Averaging	Max. 99 times
Reference value	Max. 50 pieces of reference data can be registered
External output	RS-232C
Environment for use	15 to 40°C, 30 to 80%R.H., No condensation
Power supply	100 to 240 VAC, 50/60Hz
Power consumption	100VA
Dimensions & weight	425W $\times$ 415D $\times$ 189H mm, 13kg
Related standards	JIS Z-8722, ASTM E 308, ASTM E 313, ASTM D 1925
Option	Color management software for Windows, various attachments and cells for exclusive use

※Specifications are subject to change without notice due to technical improvement.

OPTION

New Product

SE 6000 + Optical fiber, Model OF



● Spectrophotometric color difference meter



By applying a pen type head to your samples, it will become possible to measure even the ones with complicated shapes, measurement of which has been difficult before.

Specifications

Measuring diameter	Select one of 3mm ϕ , 4mm ϕ , and 6mm ϕ
Light receiving system	Select one of the following 4 systems: 0° illumination / diffused light reception, diffused illumination / 0° light reception, 0-45°, or 45-0°

Suitable for such measurements as for :

- ◇ Printed matter
- ◇ Large samples
- ◇ Convexo-concave samples
- ◇ Dental material and skin
- ◇ Automobile components
- ◇ Cosmetics
- ◇ Craftwork and timber
- ◇ Petals and leaves
- ◇ Fruits
- ◇ Beans
- ◇ Electronic components



NIPPON DENSHOKU
Advanced Technology in Color and Brightness.

⊗ NIPPON DENSHOKU INDUSTRIES CO., LTD.

HEAD SALES DEPT.

SENGOKU HASEGAWA BLDG., 4-45-17 SENGOKU, BUNKYO-KU, TOKYO 112-0011 JAPAN

PHONE : 81-3-3946-4392 - FAX : 81-3-3946-1690

URL <http://www.nippondenshoku.co.jp/>

0803-30803